


User guide for TRP-Adjustable Threshold Ramp

Before use

- Please remove all packaging and ensure the product has not been damaged in transit.
- Carefully read the instructions and warning texts before using the product
- Keep instructions for reference
- On each product you will find a label that shows the maximum permitted load.

Contents

- 1 x 14mm spanner(fig.1), 6mm Allen-key(fig.1), 2 x expansion screws (fig.2), 10 x socket head bolts (fig.2), 4 x nuts (fig.2), 4 x washers (fig.2), adjustable leg kit (fig.3), 2 leg feet (fig.3).

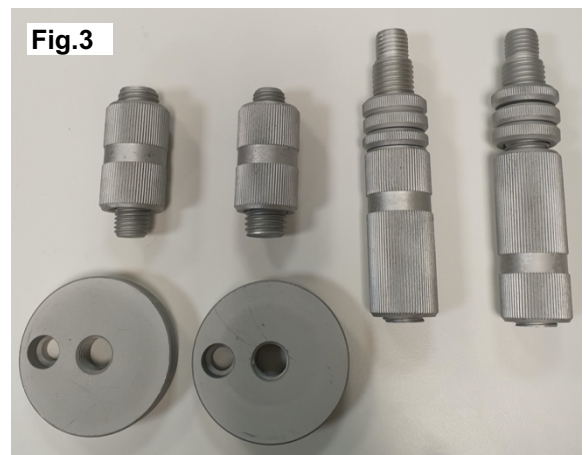

To Assemble

1. First, assemble the legs using the connections provided to your desired height. Please ensure that you use the locking nuts provided between the larger sections to ensure that they are fully secure. (fig.4)
2. Attached the ramp leg bracket to the ramp using 4 of the socket head bolts, 4 nuts and 4 washers (fig.5). Tighten using the spanner and Allen-key (fig.6). Please note that if you are using the ramp as part of the up and over kit, you will need to attach the leg bracket to the inner 4 holes on the ramp.
3. Slot the legs into the holes on the ramp let bracket (fig.7) and secure with the remaining socket head bolts provided (fig.8).

Fig. 4



Fig. 5



Fig. 6



Fig. 7



Fig. 8



Fully assembled ramp



If you wish to secure the ramp to the ground, please follow the instructions below;

1. Put the ramp in place and mark on the ground where the holes need to be drilled
2. Using an 11mm drill bit, drill holes where they are marked on the ground. Insert the bolts with the flared end with the slit at the end of the sleeve at the bottom and place the split washer over the screw and then put the ramp back in place.
3. Secure the ramp in place using the large, flat washer and nut. As this is tightened, the ramp will secure in the ground.



Maintenance

- If the product becomes damaged in any way it should not be used (contact supplier/manufacturer).
- The product may be cleaned with normal cleaning products and left to dry before use.

Warranty

The Ramp People warrants that the ramp/accessory is free from defects in material and workmanship.

This warranty shall remain in effect for 1 year from the date of original consumer purchase. Use other than the published/intended item will void the warranty.

The warranty does not cover failure of the unit during the warranty period if the failure is not caused by defect in material or workmanship or if the failure is caused by unreasonable use.

The Ramp People are not liable for incidental or consequential injury or damage to property.

Should it be necessary to return the ramp, it is essential it is adequately packed to protect it in transit, preferably in its original box; otherwise The Ramp People do not accept any responsibility for transit damage to the unit.

In some cases, the items will be supplied by 3rd party seller, please keep the original invoice. This does not affect your statutory rights.

The Adjustable Threshold Ramp is part of our range of products from The Ramp People. To view the full range please visit our website

www.theramppeople.co.uk

Or contact us via, post, email and/or telephone.

The Ramp People

Unit C3
Sandown Ind. Park Mill Road
Esher
Surrey
KT10 8BL

Tel: +(44) 1372 478960

Email: hello@theramppeople.co.uk

www.theramppeople.co.uk